IN-EXON

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
MODELLO 1.0 - 2026

ITALIANO PAG. 1-20

ENGLISH PAG. 21-37

CHINESE PAG. 38-57

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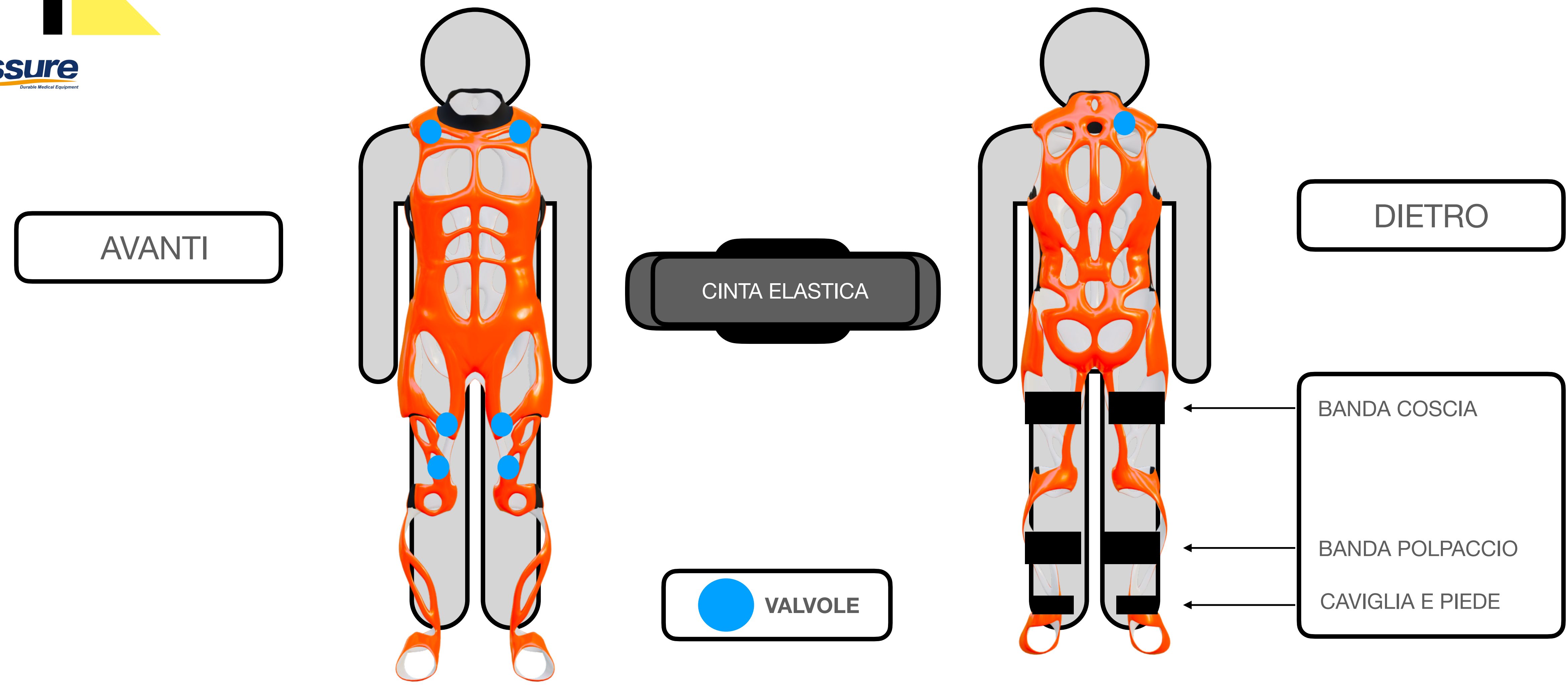- galleggiamento,
- comfort,
- performance.

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

PAG. 38-57

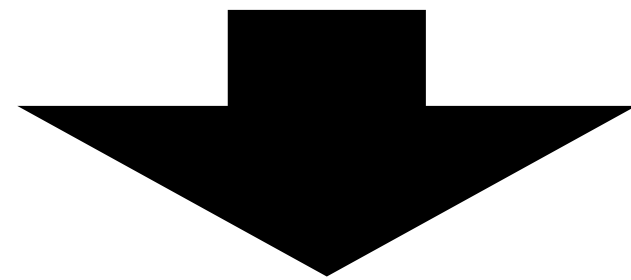
WEAR THE **IN-EXON** DEVICE

- A. DEFLATE THE DEVICE BEFORE THE USE AND AFTHET
- B. POSITION THE VELCRO OF THE INFLATABLE SHIELDS ACCORDING TO YOUR SIZE
- C. POSITION THE VELCRO OF THE LEG IN THE THIGH CONNECTION AREA
- D. CHECK ALL AIR HOSE
- E. CONNECTIONS MAKE SURE THE SHIELDS ARE PROPERLY CONNECTED


- PUT ON THE ASSEMBLED DEVICE FROM THE HEAD
- ADHESIVELY ADHERENT THE VELCRO ON THE BACK TO THE ABDOMINAL BAND
- ADHESIVELY ADHERENT THE VELCRO ON THE FRONT SHIELD
- TIGHTEN THE THIGH STRAPS
- TIGHTEN THE CALF STRAPS
- TIGHTEN THE HEEL STRAPS

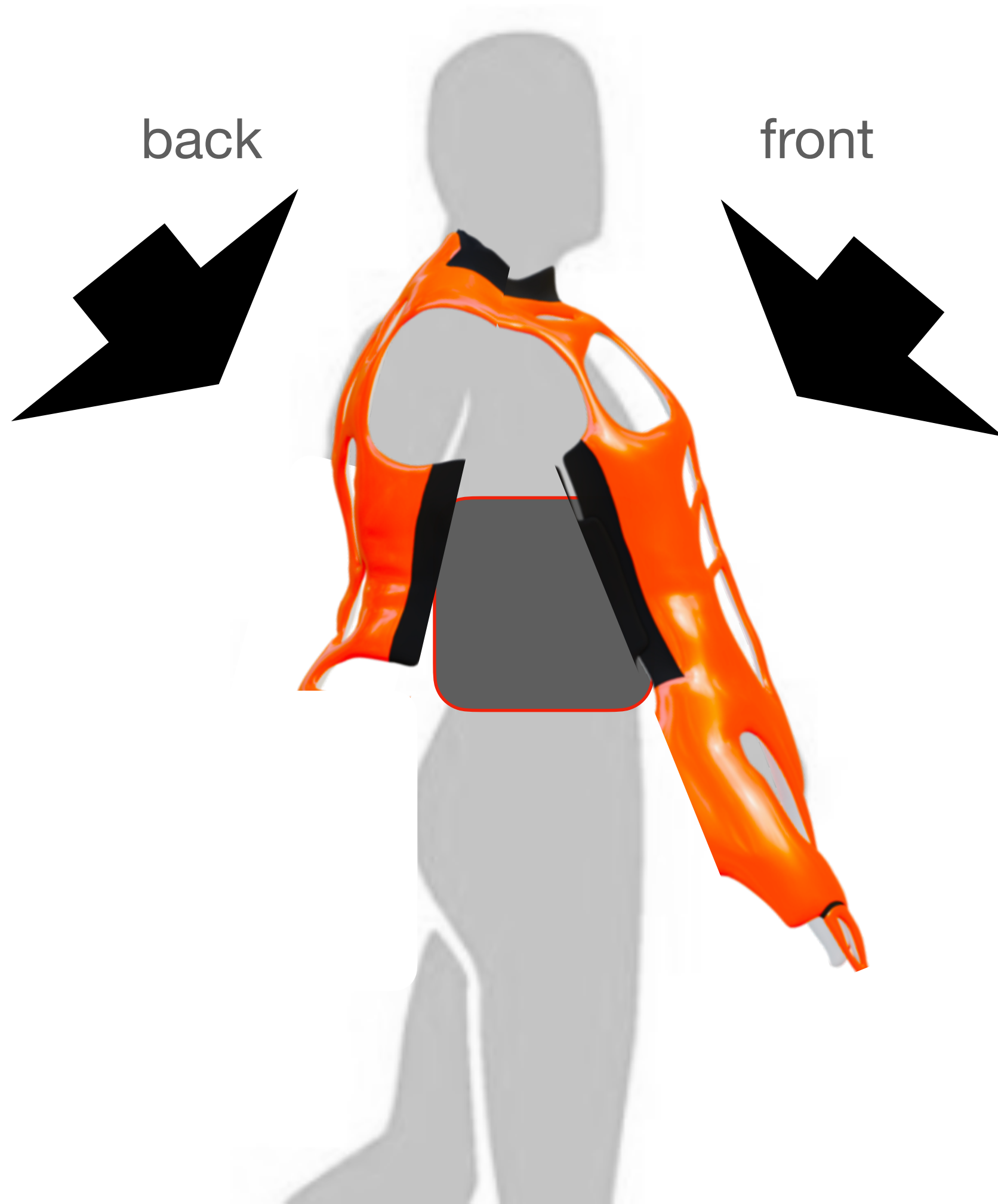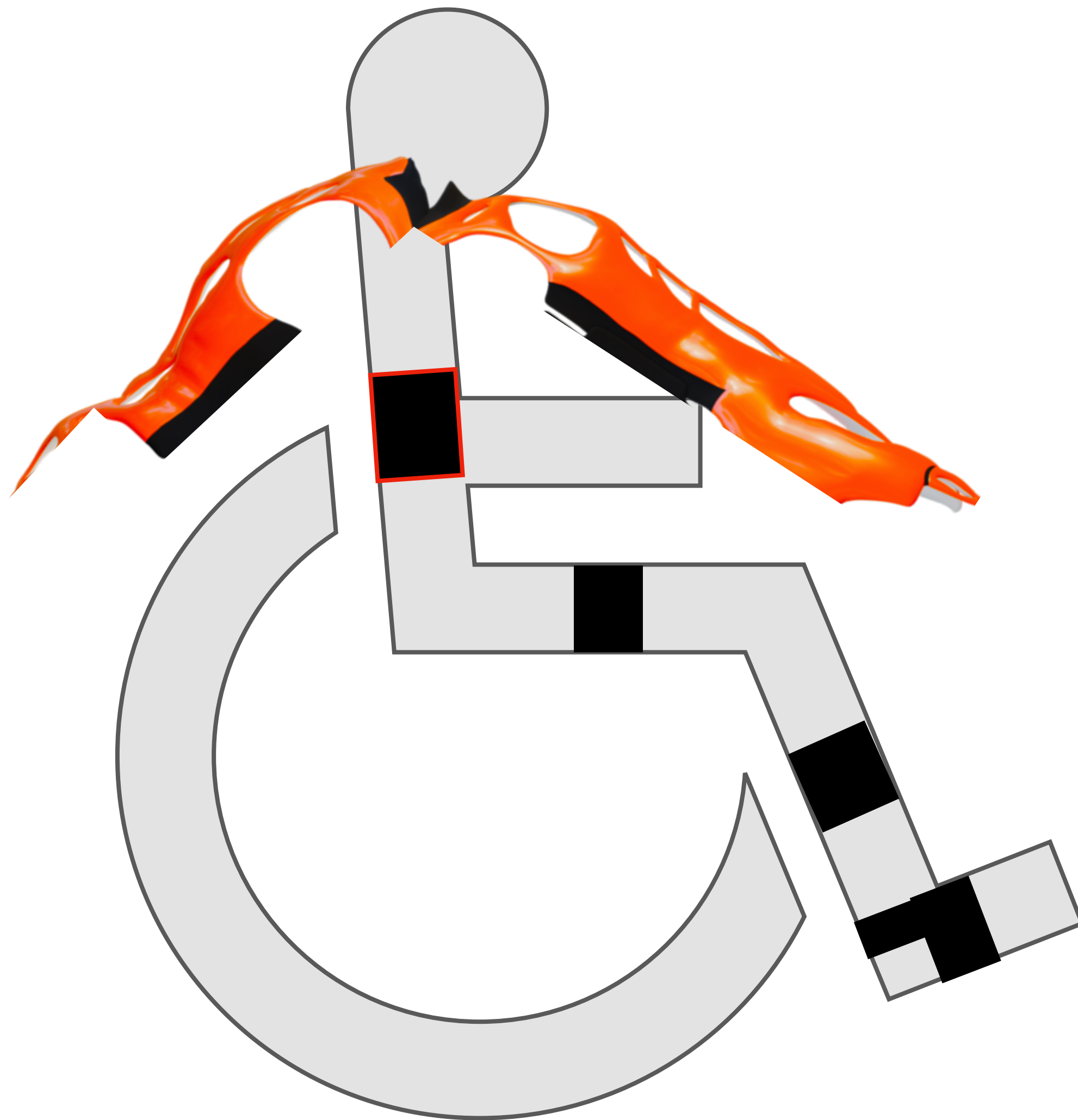IN-EXON

- 1.PUT ON THE ABDOMINAL BELT
- 2.PUT ON THE DEFLATED DEVICE
- 3.POSITION THE SHELLS AS SHOWN IN FIG. 3
- 4.TAP THE ELASTIC BANDS ON YOUR CALVES AND FEET
- 5.TAP THE ELASTIC BANDS ON YOUR THIGHS FIG. 3

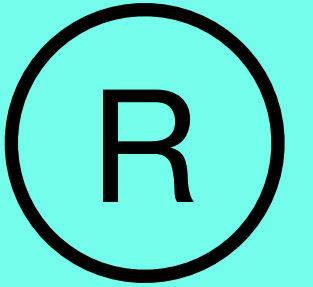ENGLISH

IN

—

EXON

Inflatable exoskeletal device
for swimming and aquatic rehabilitation

model 1.0 - 2026

FAQ

END USER

1. Is this a medical device?

No.

At present, this product is not classified as a medical device.

However, it is well suited for post-traumatic rehabilitation and assisted swimming for people with disabilities, under proper supervision.

2. What is the device used for?

It is used to:

- support buoyancy,
- improve stability in water,
- assist swimming,
- support rehabilitation and functional recovery exercises.

3. Who can use the device?

It can be used by:

- adults,
- women,
- children (only with suitable size and model),
- people with disabilities.

⚠ Children and users with reduced mobility must always be **supervised**.

FAQ

END USER

4. Where can it be used?

The device can be used:

- in swimming pools,
- in the sea,
- in supervised aquatic environments.

5. How can the device be inflated?

The device can be inflated using:

- a Bluetooth electronic pump,
- a manual pump,
- a 0.50-liter rechargeable air cylinder, using a bike-compatible accessory.

6. Can I inflate the device in the water?

Yes.

The device is designed to be:

- worn while deflated,
- inflated **directly in the water** according to personal needs.

7. How much should I inflate it?

Inflation can be adjusted to:

FAQ

END USER

8. Can I deflate it while in the water?

Yes.

The device can be **deflated in the water** to adjust:

- buoyancy,
- comfort,
- performance.

9. Can I adjust inflation multiple times during use?

Yes.

The device allows multiple inflation and deflation cycles during activity.

10. Is it safe for people with disabilities?

Yes, when used correctly.

It is suitable for **assisted swimming and water rehabilitation** but **must always be supervised**.

11. Is this device a life jacket?

No.

This device is not a life-saving device

FAQ

END USER

12. What should I do if it loses air?

Stop using the device immediately and **leave the water**.

13. How do I clean it?

- Rinse with fresh water after use.
- Use mild detergents.
- Dry completely before storage.

14. How should it be stored?

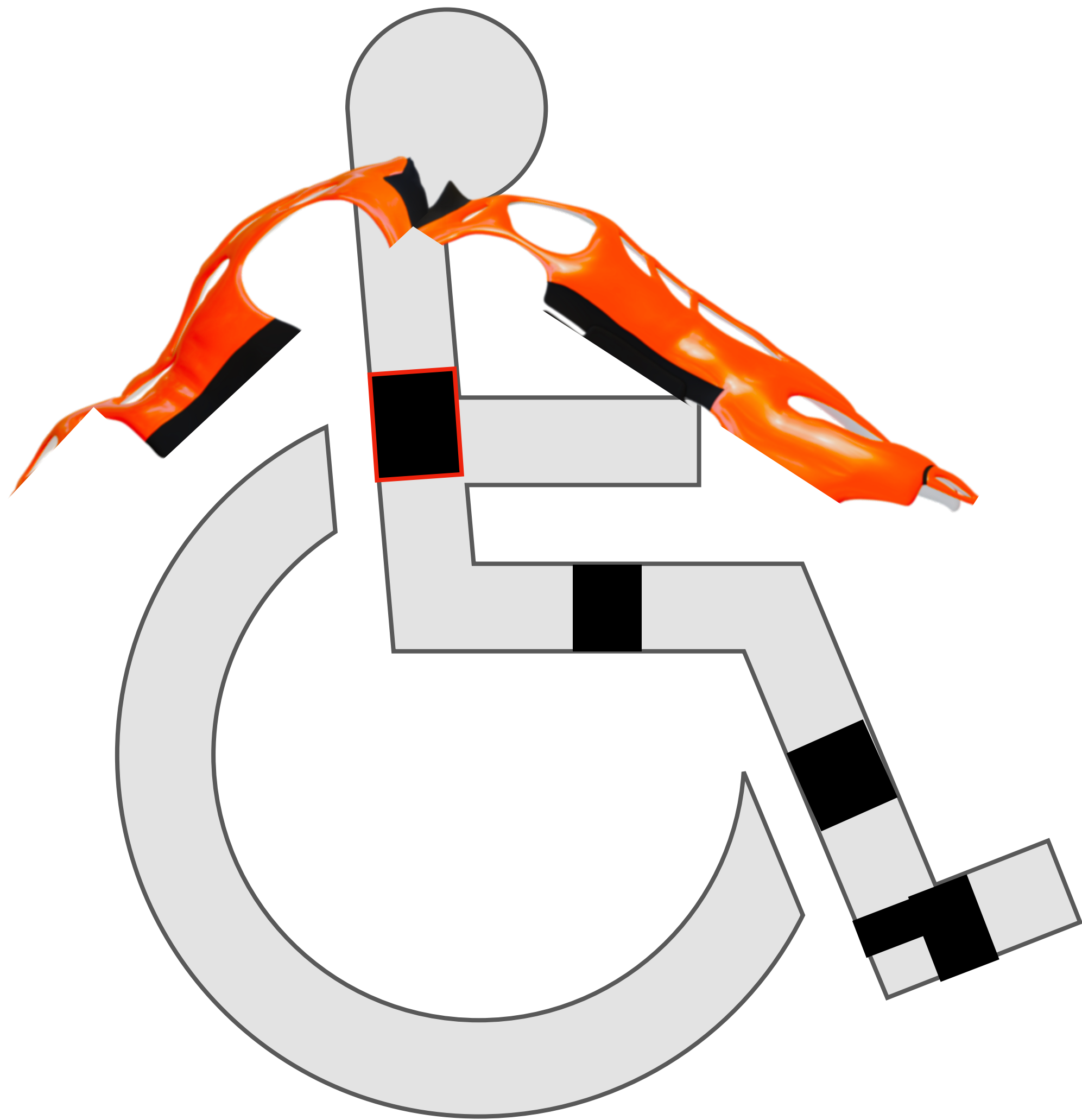
- Store **deflated**.
- In a cool, dry place.
- Away from direct sunlight and heat sources.

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